

## DESIGNED SPECIFICALLY WITH HIGH-BRIGHTNESS LEDS IN MIND: Gooch and Housego's OL 771

---





The OL 771 is a high performance, portable spectroradiometric LED test and measurement system with powerful, user-friendly software. Fast and accurate, the OL 771 is part of the most comprehensive line of LED test and measurement solutions!

## Features

The OL 771 High Speed HB-LED Measurement System is optimized to perform all critical measurements of LEDs, LED clusters, LED chips, and more. Specifically designed for fast spectroradiometric measurement of high-brightness LEDs, the system provides the necessary accuracy to characterize components for R&D purposes, while being economical and compact enough to perform rapid QC checks on the production floor.


### Features include:

- Convenient USB interface
- 25+ spectral scans/ second
- Compliant with CIE 127 guidelines with accessories for various collection geometries:
  - Averaged intensity in Conditions A & B
  - Total luminous & spectral flux
  - Goniometric profiling
- Low stray light performance
- High spectral resolution
- High dynamic range
- 0.5 nm wavelength accuracy
- Research-grade precision
- Compact, lightweight portable enclosure
- Rugged fiber optic interface with strain relief and self-centering adapter

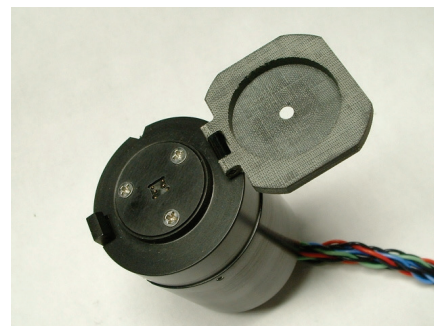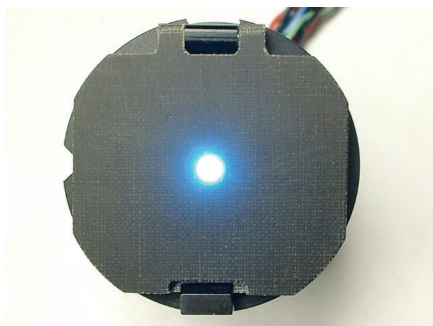

## Software

The OL 771 Application Software is highly intuitive. It shares the convenient, well-proven OL 770 software platform. The OL 771's software allows for turnkey automated operation while providing the user with easy and complete access to the system's full power and capabilities.

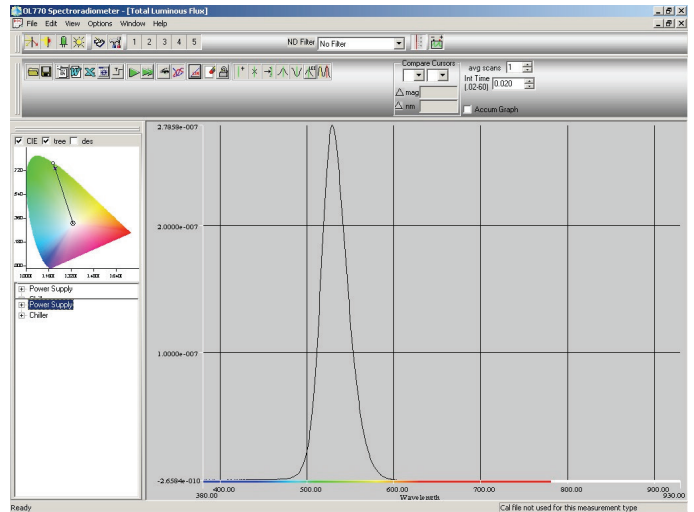


### Software features include:

- Windows XP/2000® platform
- Real-time graphics utility
- Custom report templates
- MS Excel and Word-compatible direct reporting
- Display, log, and store resultant data
- Comparison cursors
- Cursor snap to peak/ valley
- Value Monitor for real time pass/fail display of any calculated values
- Software triggering for data acquisition
- Dominant wavelength
- Peak wavelength
- Spectral bandwidth (FWHM)
- Spectral purity
- Power
- Color Temperature
- Color rendering index
- Total luminous flux
- Averaged LED intensity
- Tristimulus - 2° XYZ, 10° XYZ
- 1931 Chromaticity - 2° xyz, 10° xyz
- 1960 Chromaticity - 2° uv, 10° uv
- 1976 Chromaticity - 2° u'v', 10° u'v'
- Lab Luv - Illuminants A, B, C, D65; 2° Lab Luv, 10° Lab Luv
- Intensity profile polar plots (w/ OL 700-30)
- High-resolution CIE plots with accumulating coordinates, zooming, and dominant wavelength tracking
- Accumulation graphs for time studies
- LabVIEW example utilizing ActiveX control
- Optional ActiveX Control Software Development Kit

## Specifications

OL 771 Spectroradiometer Specifications	
Wavelength Range (standard)	380 - 780 nm
Wavelength Accuracy	± 0.5 nm
Optical Bandwidth (w/ 100 micron slit)	2.0 nm
Spectral Resolution	≈ 0.4 nm
Optical Focus Length	140 mm
Optical Input	fiber optic
Optical Aperture	f/2
Operating Temperature	0 to 30° C
Operating Humidity	0 to 90% (non-condensing)
Detector Technology	photodiode array
Integration Time	20 ms - 60 s
Dimensions	7¼" W (18.4 cm) x 13½" H (33.6 cm) x 13" D (33.0 cm)
Weight	21.0 lbs (9.5 kg)

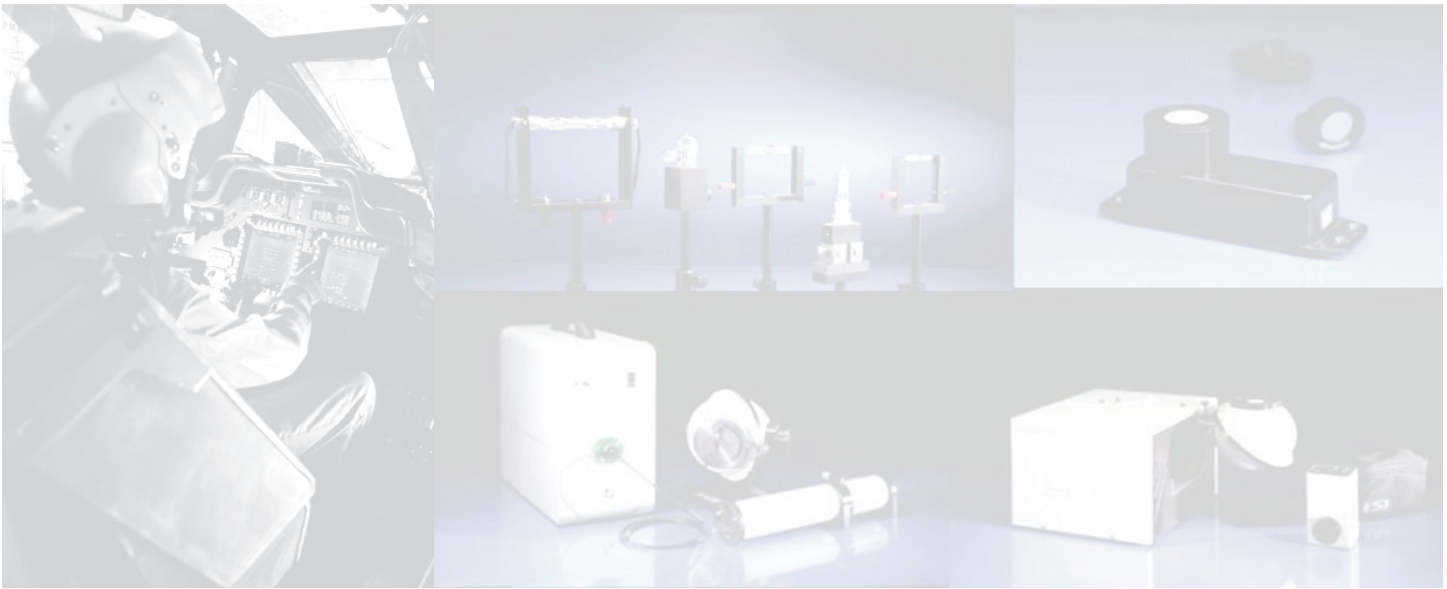

All Gooch & Housego products are supported by outstanding factory customer service and a worldwide network of representatives and distributors. Put our 40 years of expertise in precision light measurement to work for you.

Night vision cockpit images used with permission from United Rotocraft.